

AUTOSAMPLERS & SAMPLE COLLECTION

CD AUTOPLUS™ & AUTOFILL™




The CD AutoPlus™ is a precision autosampling system with all-inert syringes, valves, tubing, and fluidics, designed for plug-and-play integration with CD14 dissolution testers. With 6 syringes (10 mL capacity) and a volume accuracy of $\pm 1\%$ (0.1 mL), the CD AutoPlus easily accommodates sampling intervals under 5 minutes.

Use the DissoScan™ configuration with optional return-to-vessel (RTV) probes to recirculate rinse volumes, or the Maximizer™ for automated addition of dilutions or stabilizers, and its included media-replacement kit with low-evaporation vessel covers, 6 fixed probes (PEEK), and a 6-channel tubing harness.

The CD AutoPlus (38 cm W x 66 cm D) when combined with the CD14 dissolution tester provides the smallest 14-vessel footprint in its class.


- Handles surfactants and buffered media
- Sample timing options: At Point or Midpoint
- Seamless integration with space-saving AutoFill™ fraction collector stacked on top
- Convenient auto-start feature
- Programmable end-of-test wash routine automatically prepares system for next test
- Print to serial and network printers
- 21 CFR Part 11 security compliant with event logs, error logs, and full audit trail
- USP / US FDA / ASTM / EP / JP / CE / CSA / RoHS compliant

CD AUTOPLUS™ & AUTOFILL™

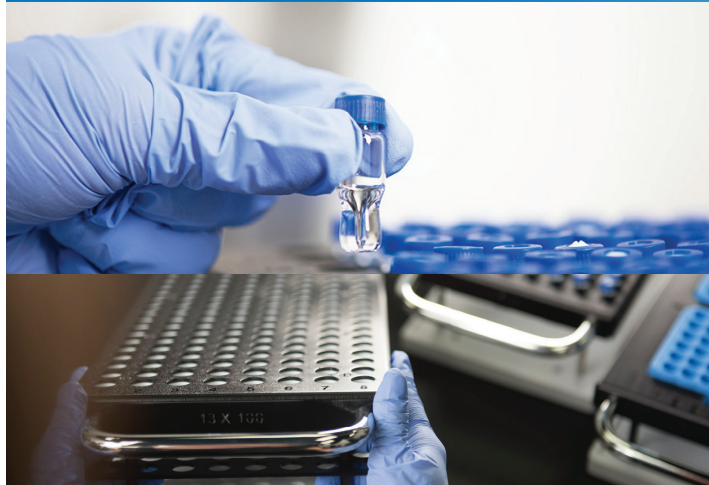


AUTOFILL™ FEATURES

- Precision fraction collector fully compatible with all AutoPlus™ autosampler models
- Wash port with end-of-test wash routine
- Sample retrieval return to vessel accessory
- CE / CSA / RoHS compliant
- Choice of easy install collection racks with rack detection sensors
- Compatible with 16 x 100 mm test tubes with evaporation caps, and 12 x 32 mm vials with pre-cut septa

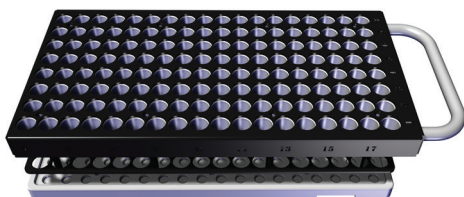

Collection Rack for 16 x 100 mm Test Tubes

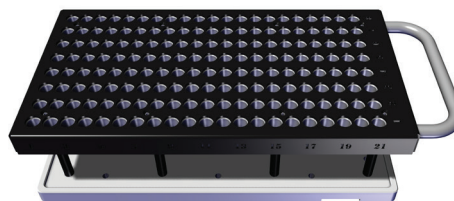

Collection Rack for 12 x 32 mm HPLC Vials

ORDERING INFO

For more information, please contact your local Teledyne Hanson representative or email us at: hansonsales@teledyne.com

CD AUTOPLUS MAXIMIZER

61-101-206 CD AutoPlus Maximizer, 6-Sample, 115/230 V

AUTOFILL FRACTION COLLECTOR

61-201-001 AutoFill Collector, 115/230 V

CD AUTOPLUS DISSOSCAN

61-101-216 CD AutoPlus DissoScan, 6-Sample, 115/230 V

COLLECTION RACKS

61-201-052 Collection Rack for 16 x 100 mm Test Tubes (18 rows)

61-201-053 Collection Rack for 12 x 32 mm HPLC Vials (21 rows)



9810 Variel Avenue
Chatsworth, CA 91311 USA

+1 818.882.7266
www.teledynehanson.com

© Teledyne Hanson, a division of Teledyne Instruments, Inc. | 99-310-075 Rev. 26 Oct 2020