


AUTOSAMPLERS & SAMPLE COLLECTION

VISION[®] AUTOPLUS[™] & AUTOFILL[™]

The Vision[®] AutoPlus[™]

is a modern upgrade of the popular Hanson legacy AutoPlus[™] autosampler. The Vision AutoPlus comes with a choice of 6 or 8 sampling channels, and DissoScan[™] (1-bath) or Maximizer[™] (up to 3 baths) configurations. Precision syringe pump mechanics, all-inert tubing and valves, and color touchscreen programming with up to 100 protocols are standard on each instrument.

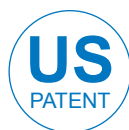

The Vision[®] AutoFill[™]

is a sample collector designed for seamless integration with the Vision AutoPlus. Vision AutoFill accessories include choice of collection racks with 13 x 100 or 16 x 100 mm test tubes for UV-Vis, 12 x 32 mm vials for HPLC, sample trays with 12 x 32 mm vials for UPLC[®], and more.

VISION® AUTOPLUS™ & AUTOFILL™



VISION® AUTOPLUS™ FEATURES

- Precision syringe pump autosampling system
- All-inert syringes, valves, tubing, and fluidics
- Choice of 6 or 8 sampling syringes (for 6 or 8 vessels)
- Volume accuracy and repeatability $\pm 1\%$
- Seamless integration with G2 Classic 6™ and Elite 8™
- Automatic media replace standard with Maximizer™
- Convenient “Auto Start” feature
- Automatic addition of dilutions or stabilizers
- Optional in-line filter block kit (easy tubing and filter connections)
- Handles problem drugs, surfactants, and buffers
- Trouble-free maintenance, no peri-pumps or check valves
- Small, stacked footprint leaves open bench space in front for waste beaker and other tools
- Customizable end of test auto-wash feature automatically prepares the system for the next test
- Full color touchscreen with intuitive menus and programming
- Up to 100 protocols (create, manage, and save routines)
- Upload and download programs with flash drive
- Electronics and firmware include 32-bit digital technology
- Print to both serial and network (PostScript) printers
- Alarms for sampling, maintenance, and calibration
- Compatible with non-Hanson dissolution testers
- Security: 50 users
- 21 CFR Part 11 security compliant
- USP / US FDA / ASTM / EP / JP / CE / CSA / RoHS compliant

VISION® AUTOFILL™ FEATURES

- Precision fraction collector built for the Vision® platform
- Seamless integration with AutoPlus™
- Choice of easy install collection racks with rack detection sensors
- Collect in test tubes for UV or vials for HPLC & UPLC®
- Compatible with 13 x 100 and 16 x 100 mm test tubes and 12 x 32 mm vials
- New wash port with automatic end-of-test wash routine
- Sample retrieval return to vessel accessory
- CE / CSA / RoHS compliant

ORDERING INFO

For more information, please contact your local Teledyne Hanson Research representative or email us at: hansonsales@teledyne.com

- 61-101-006** Maximizer, 6-Sample, 115/230 V
- 61-101-008** Maximizer, 8-Sample, 115/230 V
- 61-101-016** DissoScan, 6-Sample, 115/230 V
- 61-101-018** DissoScan, 8-Sample, 115/230 V
- 61-201-001** AutoFill Collector, 115/230 V